


# User Manual for Automated Lever Task Arena

designed by Aiva levins  
script created by Adrien Boissenin  
published by Soshi Samejima

soshis@uw.edu

copyright 2021

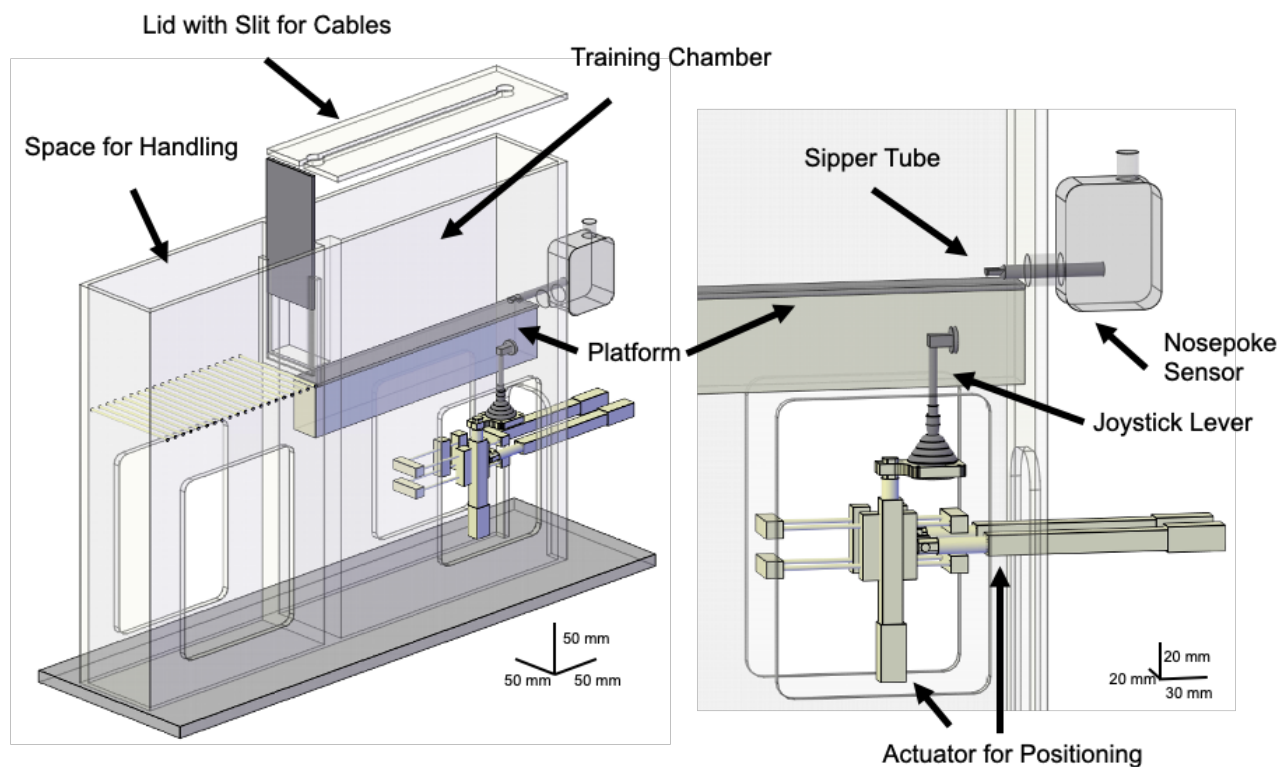

## File List:

Materials: List of needed materials and dimensions

Published paper

Supplemental Video

## Folder List:

code : python code to launch desktop logging program and run steps

build-files: \*.dwg files for 3D printed items and the arena design

## Use training python script:

Open 'Training\_Step(number).py' file in lever arena folder -> press F5 (should make the program start) -> type in animal's number -> start training -> After training press Ctrl+C to terminate session.

Exception for Step6: Open "Training\_Step(number).py" file in lever arena folder -> press F5 (should make the program start) -> type in animal's number -> type last vertical position reached (0, 1 or 2) -> type last horizontal position reached (0, 1, 2, 3 or 4) -> start training -> After training press Ctrl+C to terminate session.

## \*\*\* Phase 1 : nose spoke, lever push.

- Step 1: perform nose poke. Success for nose poke is notified with beep sound and reward. After a single session with 40 trials, the step is completed. Don't forget to progressively adjust the gain knob to its minimum. Nose spoke should not be neglected. Rats that don't learn to do nose spoke right inside the sensor neck performed poorer when the task become more complicated. They learn slower as they get easier confuse and their lever pushes are less 'proper'(expected time: 2 sessions).

- Step 2: Learning to push the lever. Each time, start with 5 nose spokes to get the reward. First real lever is in low position ( $v\_pos = 0.00110$ ,  $h\_pos = 0.0130$ ). Use dummy lever to draw rat's attention. Give reward when rat touches the dummy lever and then when the rat touches the dummy lever. Start by putting the dummy lever at a similar level than the platform and then, to prevent him

using its teeth, lower your dummy lever. Once rat learned to push a little bit the dummy lever, shift attention from dummy lever to real lever by rising the real lever to  $h\_pos = 0.00120$  and finally  $h\_pos = 0.00130$  (threshold = 2500). The script will stop and ask if you want to rise the lever after the 10th, 20th and 30th successes. press 'y' (yes) or 'n' (no) then press 'enter' to move it -or not-. Once you believe rats acquire the movement step is complete. (Expected time: 2 sessions).

- Step 3: Learning to push the lever on a wider range with an adaptive threshold. Joystick is moved at  $v\_pos = 0.00130$  and  $h\_pos = 0.00120$  horizontal position (threshold set at 2400). Success for the trial is notified with beep sound. Lever threshold increases each 5 corrects push. If 5 incorrect pushes on the Lever before a successful trial, the lever threshold decrease. When the rat reached a threshold of at least 3100 (maximum threshold is 3200), the step is complete. (Expected time: 1 sessions or 2 sessions, can be done right after step2 if training time remains).

- Step 4: Association nose spoke and lever push: Begin with 10 nose spokes-only with a lever at an unreachable position ( $v\_pos = 0.00110$ ). Then, lever is moved:  $v\_pos = 0.00130$  and  $h\_pos = 0.00120$ , threshold= 2400. Usually, the rat will push the lever, maintaining it pushed while exploring the environment with its nose. once he learns the association: push + nose spoke => reward, he will become more efficient (doing nose spoke, then a short push on the lever). Again, don't underestimate the importance of the step. The rat needs to stick its head right in front of the sensor. Rat performing well here, usually learn quicker the most complex step as they can focus on the lever push. When you noticed the rat learn to do both tasks to get the reward, go step 5. (Expected time: 2 sessions).

\*\*\*Phase 2: Maximum lever amplitude at every position of the lever.

- Step 5: Last and longest step. In this step, the lever is set at the initial position ( $v\_pos = 0.00130$  and  $h\_pos = 0.00105$ ) and will move to every position row after row. The rat is expected to push the lever, beyond the threshold of 3200, 4 times before moving to the next position. An adaptive threshold is setup to assist the rat in that task. First, 4 incorrect pushes allow to reduce the threshold, the number is low to keep the rats motivated to do a large number of pushes. However, for the last row, it was observed that the rats performed better when the number of pushes allowing to reduce the threshold is set to 7 or 8. The task is considered complete when the rats reach maximum threshold (3200) 4 times at row 2 and column 4 (most remote position). (Expected time: as long as needed, 4-8 sessions).

- Step 6: Last step. Starting from row 0, column 0, the goal is to go through each position 2 times (aka 90 successes). After three corrects push the lever will move

(threshold at 3200). Once complete the training is over. I observed that to be completed the first time the rat needs at least 20 minutes but this time is then reduced to 15 min the next session (expected time: 1 or 2 sessions).

\*\*\*Phase 3: After the training

- Step Random: Once the animals are trained, you can reduce the training to 2 or 3 5min-session per week. At this step, the lever moves to a random position (among the 15 learned) every 3 correct pushes). Usually maintain the training until you stop the water restriction before the surgery.